**SUNRISE
MEDICAL**

EASY 200/300 **EASY** MAX.

GEBRAUCHSANLEITUNG
DIRECTIONS FOR USE
NOTICE D'UTILISATION
MANUALE PER L'USO
GEBRUIKSHANDLEIDING
MANUAL DE USUARIO

Vorwort (1.1.1)

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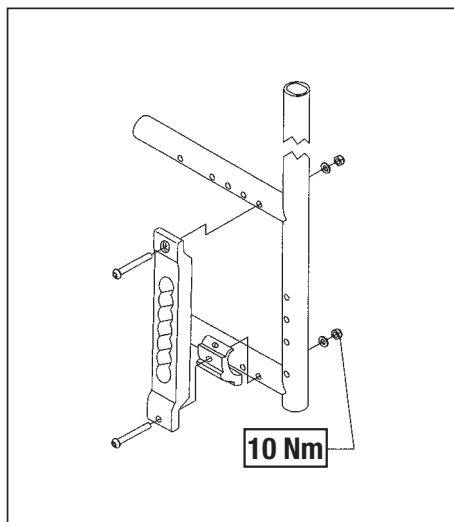
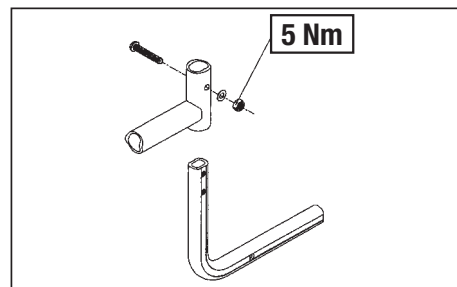
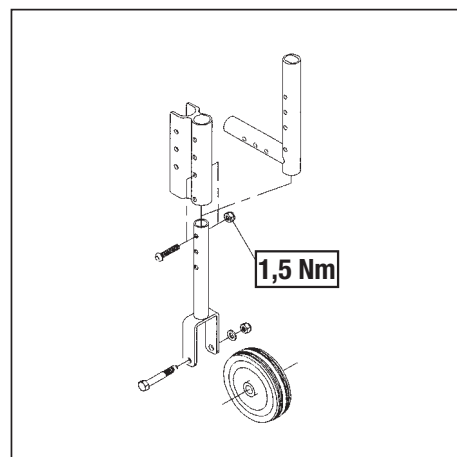
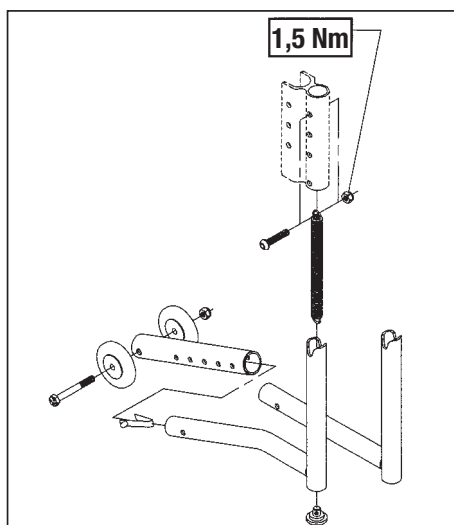
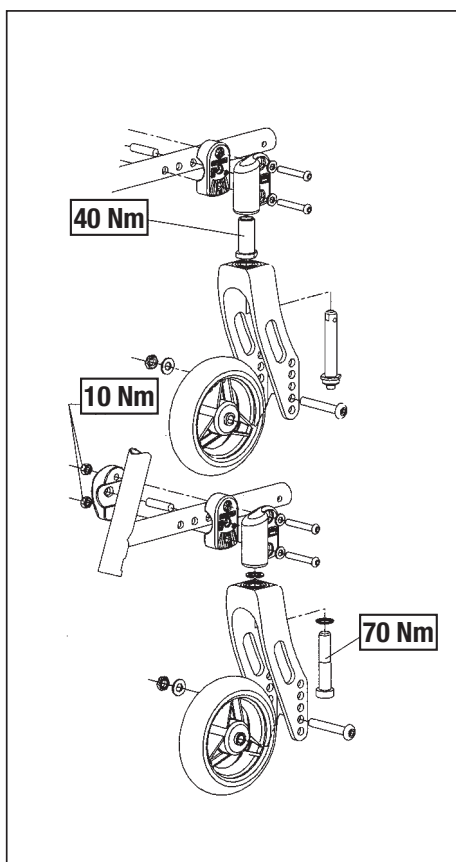
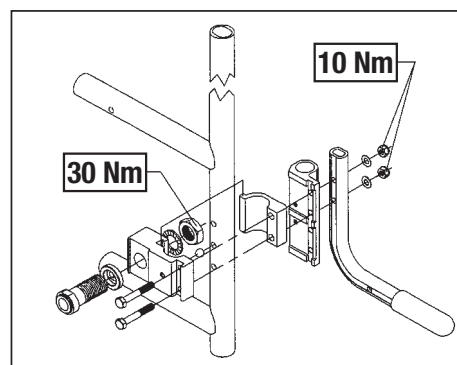
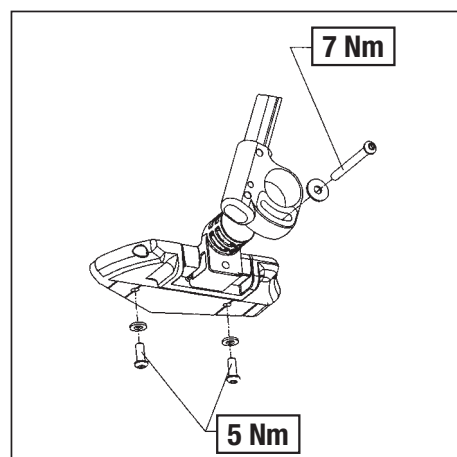
Foreword (1.1.1)

Dear Customer,

We are very happy that you have decided in favour of a high-quality product from SUNRISE MEDICAL.

This user's manual will provide numerous tips and ideas so that your new wheelchair can become a trustworthy and reliable partner in your life.

Maintaining close links with our customers is of great importance to us at Sunrise Medical. We therefore like to keep you up-to-date with our new and current developments. Keeping **close** to our customers also means fast service when you need replacement parts or accessories, or just have a question about your wheelchair – and with as little red tape as possible.

We want you to be satisfied with our products and service. Sunrise Medical therefore constantly works at continuous development of its products. For this reason, changes can occur in our range of products with regard to shape, technology, and fittings. Consequently, no claims can be construed from the data or pictures contained in this user's manual.

SUNRISE MEDICAL has been awarded the ISO 9001 Certificate, which affirms the quality of our products at every stage, from R & D to production.

Please contact your local, authorised SUNRISE MEDICAL dealer if you have any questions concerning the use, maintenance, or safety of your wheelchair.

As the manufacturer, SUNRISE MEDICAL, declares that the lightweight wheelchairs conform to the 93/42/EEC / 2007/47/EEC guideline.

In the case that there is no authorised dealer in your area or you have any questions, you can contact Sunrise Medical either in writing or by telephone at:

**Sunrise Medical Ltd.
Sunrise Business Park
High Street, Wollaston
West Midlands DY8 4PS
England**

Telephone: +44/1384-446688

Fax: +44/1384-446699

www.sunrisemedical.co.uk

IMPORTANT:

DO NOT USE YOUR WHEELCHAIR UNTIL THIS MANUAL HAS BEEN READ AND UNDERSTOOD.

Table of contents

Foreword for Wheelchairs

General safety notes and driving restrictions18-19

Wheelchair Components19

Handling

Folding Up and Unfolding20

Options

Tipping Tubes20

Brakes20-21

Footplates21

Footrest22

Headrest22

Casters22

Axle plates23

Cambers23

Amputee Axle Plate23

Backrests23

Sideguards23-24

Push Handles24

Lap belt instructions24

Anti-tip Tubes25

Seat Sling25

Crutch Holder25

Tray25

Stabilizing Bar25

Transit Wheels26

Tires and Mounting26

Possible Problems26

Maintenance26

Nameplates/Guarantee26

Technical Data27-28

Torque29

Use

Wheelchairs are exclusively for a user who is unable to walk or has limited mobility, for their own personal use in- and outdoor.

The maximum weight limit (includes the user and any weight of accessories fitted to the wheelchair) is marked on the serial number label, which is affixed to the crossbar or stabiliser bar below the seat.

Warranty can only be taken on if the product is used under the specified conditions and for the intended purposes.

The expected life of the wheelchair is 5 years.

Please **DO NOT** use or fit any 3rd party components to the wheelchair unless they are officially approved by Sunrise Medical.

Area of application

The varieties of fitting variants, as well as the modular design, mean that it can be used by those who cannot walk or have limited mobility because of:

- Paralysis
- Loss of extremity (leg amputation)
- Extremity defect/deformity
- Joint contractures/joint injuries
- Illnesses such as heart and circulation deficiencies, disturbance of equilibrium or cachexia as well as for elderly people who still have strength in the upper body.

When considering provision, please also note the body size, weight, physical and psychological constitution, the age of the person, living conditions and environment.

General safety notes and driving restrictions (3.2)

The engineering and construction of this wheelchair has been designed to provide maximum safety. International safety standards currently in force have either been fulfilled or exceeded in parts. Nevertheless, users may put themselves at risk by improperly using their wheelchairs. For your own safety, the following rules must be strictly observed.

Unprofessional or erroneous changes or adjustments increase the risk of accident. As a wheelchair user, you are also part of the daily traffic on streets and pavements, just like anyone else. We would like to remind you that you are therefore also subject to any and all traffic laws.

Be careful during your first ride in this wheelchair. Get to know your wheelchair.

Before each use, the following should be checked:

- Quick-release axles on the rear wheels
- Velcro on seats and backrests
- Tyres, tyre pressure and wheel locks.

Before changing any of the adjustments of this wheelchair, it is important to read the corresponding section of the user's manual.

It is possible that potholes or uneven ground could cause this wheelchair to tip over, especially when riding uphill or downhill.

When riding over a step or up an incline frontally, the body should be leaning forward.

Danger!

NEVER exceed the maximum load of 125 kg for driver plus any items carried on the wheelchair. If you exceed the maximum load, this can lead to damage to the chair, or you may fall or tip over, lose control and may lead to serious injury of the user and other people.

Danger!

When it is dark, please wear light clothing or clothing with reflectors, so that you can be seen more easily. Make sure that the reflectors on the side and back of the wheelchair are clearly visible. We would also recommend that you fit an active light.

Danger!

To avoid falls and dangerous situations, you should first practice using your new wheelchair on level ground with good visibility.

Danger!

When getting on or off the wheelchair, do not use the footboards. These should be flipped up beforehand and swung to the outside as far as possible. Always position yourself as close as possible to the place where you wish to sit.

Danger!

Only use your wheelchair properly. For example, avoid travelling against an obstacle without braking (step, kerb edge) or dropping down gaps.

Danger!

The wheel locks are not intended to brake your wheelchair. They are only there to ensure that your wheelchair does not begin rolling unintentionally. When you stop on uneven ground, you should always use the wheel locks to prevent such rolling. Always apply both wheel locks; otherwise, your wheelchair could tip over.

Danger!

Explore the effects of changing the centre of gravity on the behaviour of the wheelchair, for example on inclines, slopes, all gradients or when overcoming obstacles. Do this with the secure aid of a helper.

Danger!

With extreme settings (e.g. rear wheels in the most forward position) and less than perfect posture, the wheelchair may tip over even on a level surface.

Danger!

Lean your upper body further forward when going up slopes and steps.

Danger!


Lean your upper body further back when going down slopes and steps. Never try to climb and descend a slope diagonally.

Danger!

Avoid using an escalator which may lead to serious injury in the event of a fall.

Danger!

Do not use the wheelchair on slopes $> 10^\circ$. The Dynamic safe slope is dependant on the chair configuration, the users abilities and the style of riding. As the users abilities and style of riding cannot be pre-determined then the max safe slope cannot be determined. Therefore this must be determined by the user with the assistance of an attendant to prevent tipping. It is strongly recommended that inexperienced users have Ant-tips fitted.

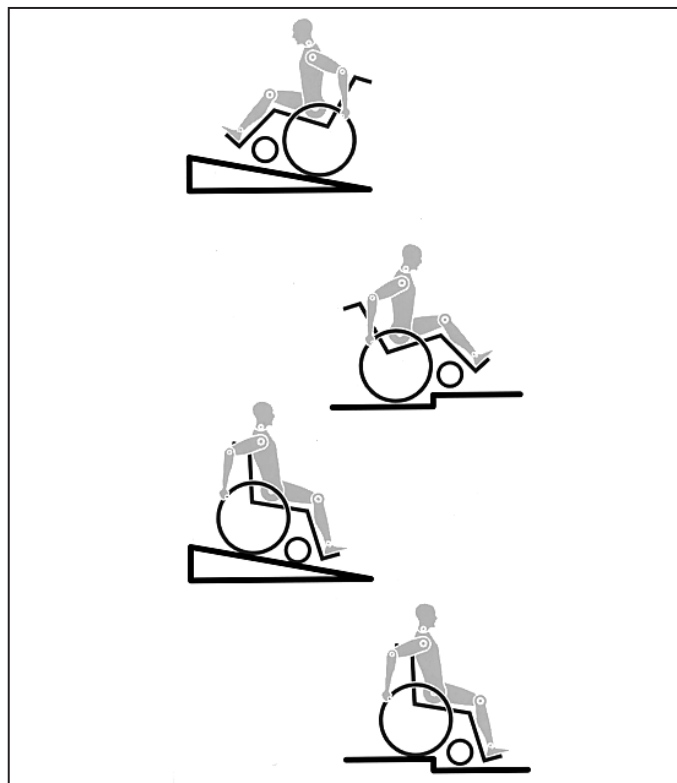Danger!

It is possible that potholes or uneven ground could cause this wheelchair to tip over, especially when riding uphill or downhill.

Danger!

Do not use your wheelchair on muddy or icy ground. Do not use your wheelchair where pedestrians are not allowed.

Danger!

To avoid hand injuries do not grab in between the spokes or between the rear wheel and wheel lock when driving the wheelchair.

Danger!

In particular when using lightweight metal handrims, fingers will easily become hot when braking from a high speed or on long inclines.

Danger!

Only attempt stairs with the help of an attendant. There is equipment available to help you, e.g. climbing ramps or lifts, please use them. If there is no such equipment available, then the wheelchair must be tipped and pushed, never carried, over the steps (2 helpers). We recommend that users over 100 kg in weight do not use this stairway manoeuvre!

In general, any anti-tip tubes fitted must be set beforehand, so that they cannot touch the steps, as otherwise this could lead to a serious tumble. Afterwards the anti-tip tubes must be set back to their correct position.

Danger!

Make sure that the attendant only holds the wheelchair using securely mounted parts (e.g. not on the footrests or the sideguards).

Danger!

When using the lifting ramp make sure that the anti-tip tubes fitted are positioned outside the danger area.

Danger!

Secure your wheelchair on uneven ground or when transferring (e.g. into a car) by using the brakes.

Danger!

If and whenever possible, during a journey in a specially fitted vehicle for disabled people, vehicle occupants should use the seats in the vehicle and the appropriate restraint system. This is the only way to ensure that occupants will have the maximum protection if there is an accident. When using safety elements offered by SUNRISE MEDICAL and using a specially designed safety system, lightweight wheelchairs can be used as a seat when being transported in a specially fitted vehicle. (See the Chapter on "Transportation").

Danger!

Depending on the diameter and setting of the castors, as well as

the centre of gravity setting of the wheelchair, the castors may begin to wobble at high speeds. This can lead to the castors being blocked and the wheelchair may tip over. Therefore, please make sure that the castors are adjusted correctly (see the Chapter "Castors"). In particular, do not travel on an incline without brakes, travel at a reduced speed. We recommend that novice users use anti-tip tubes.

Danger!

Anti-tip tubes should prevent the chair tipping over backwards unintentionally. Under no circumstances should they take the place of transit wheels, and be used to transport a person in a wheelchair with the rear wheels removed.

Danger!

When reaching for objects (which are in front of, to the side or behind the wheelchair) make sure that you do not lean too far out of the wheelchair, as if you change the centre of gravity there is a risk of tipping or rolling over. The hanging of additional load (back pack or similar items) onto your chair backposts can affect the rearward stability of your chair, especially when used in combination with recliner backrests. This can cause the chair to tip backwards causing injury.

Danger!

For thigh amputees you must use anti-tip tubes.

Danger!

Before setting off, check that your tyre pressure is correct. For rear wheels it should be at least 3.5 bar (350 kPa). The max. pressure is indicated on the tyre. The knee-lever brakes will only work if there is sufficient tyre pressure and if the correct setting has been made (see the Chapter on "Brakes").

Danger!

If the seat and back sling are damaged, you must replace them immediately.

Danger!

Be careful with fire, in particular with burning cigarettes. Seat and back slings can be set alight.

Danger!

If the wheelchair is subject to direct sunlight for a long period of time, then parts of the wheelchair (e.g. frame, legrests, brakes and sideguard) may become hot (>41°C).

Danger!

Always make sure that the quick-release axles on the rear wheels are set properly and lock in. If the button on the quick-release axle is not pressed in, the rear wheel cannot be removed.

Note !

The effect of the knee-lever brake as well as the general driving characteristics are dependant on tyre pressure. The wheelchair is significantly lighter and easier to manoeuvre when the rear wheels are pumped up correctly and both wheels have the same pressure.

Note !

Make sure that your tyres have sufficient tread! Please note that you are subject to any and all traffic laws when driving in public traffic.

Note !

Always be careful with your fingers when working or adjusting the wheelchair!

The products shown and described in this manual may not be exactly the same in every detail as your own model. However, all instructions are still entirely relevant, irrespective of detail differences. The manufacturer reserves the right to alter without notice any weights, measurements or other technical data shown in this manual. All figures, measurements and capacities shown in this manual are approximate and do not constitute specifications.

Transportation

Your wheelchair is not suitable for use as a seat in vehicles.

Do not sit in the wheelchair during transportation. Use the vehicle seat or a properly secured transport approved seat.

Wheelchair Components (4.1)

We at SUNRISE MEDICAL have been awarded the ISO 9001 Certificate, which affirms the quality of our products at every stage, from R & D to production. This product complies with the standards set forth in EU directives. Optional equipment and accessories are available at extra charge.

Wheelchairs:

1. Push handles
2. Backrest sling
3. Sideguards
4. Seat sling
5. Footrest
6. Casters
7. Caster plates
8. Fork
9. Quick-release axles
10. Wheel locks
11. Handrim
12. Rear wheel

Handling

Folding up (6.1)

First remove the seat cushion from the wheelchair and flip up the (platform) footplate or the individual footplates. Take hold of the sling from underneath in the middle and pull upwards. The wheelchair will then fold up. To fold your wheelchair until it is as small as possible, e.g., to put into a car, you can remove the footrests (depending on the model). For this purpose, open the latch from the outside and swing the footrest to the side. Then the footrest only has to be drawn out of the frame tubing.

02.1

02.2


Unfolding

Open the strap on folded frame prior to unfolding.

Press down on the seat tubing (see picture). Your wheelchair is then unfolded. Now snap the seat tubing into position in the seat saddle. This can be facilitated by slightly tipping your wheelchair. In this way, the weight can be taken off of one rear wheel so that the wheelchair can be more easily unfolded. Be careful not to get your fingers caught in the cross-tube assembly.

Quick-Release Axles for Rear Wheels (6.2)

The rear wheels are equipped with quick-release axles. The wheels can thus be installed or removed without using tools. To remove a wheel, simply depress the quick-release button on the axle (1) and pull it out.

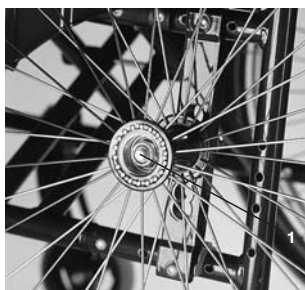04

CAUTION:

Hold the quick-release button on the axle depressed when inserting the axle into the frame to mount the rear wheels. Release the button to lock the wheel in place. The quick-release button should snap back to its original position.

Quick-Release Axles for Casters (6.3)

Castors can also be equipped with quick-release axles (1), which can then either be installed or removed without using tools.

06

Options – Tipping Tubes

Tipping Tubes (7.1)

Tipping tubes are used by attendants to tip a wheelchair over an obstacle. Simply step on the tube to push a wheelchair, for example, over a curb or step.

07

Options – Brakes

Wheel Locks (7.10)

Your wheelchair is equipped with two wheel locks. They are applied directly against the tires. To engage, press both wheel-lock levers forward against the stops (Fig. 10.1 and 10.3). To release the wheel, pull the levers back to their original positions.

Braking power will decrease with:

- Worn tire tread
- Tire pressure that is too low
- Wet tires
- Improperly adjusted wheel locks.


The wheel locks have not been designed to be used as brakes for a moving wheelchair. The wheel locks should therefore never be used to brake a moving wheelchair. Always use the hand-rims for braking. Make sure that the interval between the tires and wheel locks complies with given specifications (Fig. 10.2 and 10.4). To readjust, loosen screw (1) and set the appropriate interval. Tighten screw (see the page on torque)

CAUTION:

After each adjustment of the rear wheels, check the interval to the wheel locks and readjust if necessary.


10.1 Knee-lever brake

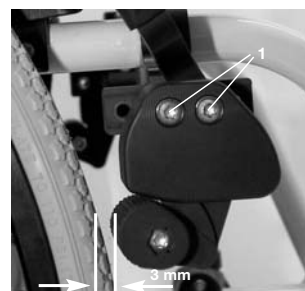10.2 Knee-lever brake

10.3 Standard brake

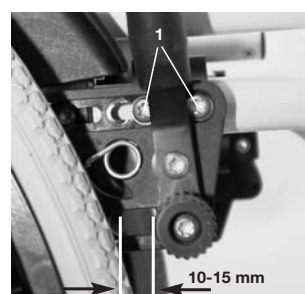10.4 Standard brake

Options – Brakes

Extension for the Brake Lever (7.11)

The extension for the brake lever can be removed or folded down. The longer lever helps to minimize the effort needed to set the wheel locks.

Caution:

Mounting the wheel lock too close toward the wheel will result in a higher effort to operate. This might cause the brake lever extension to break!

Leaning onto the brake lever extension while transferring will cause the lever to break! Splashing water from tires might cause the wheel lock to malfunction.

11

Drum Brakes (7.13)

Drum brakes permit safe and convenient braking for an attendant. They can also be set with the aid of a locking lever (1) to prevent rolling. This lever must snap audibly into place. Drum brakes do not independent on the aid of air pressure inside the tires. Your wheelchair cannot be moved when the drum brakes are set.

12

CAUTION:

Drum brakes should only be adjusted by authorized dealers.


One Hand Wheel Locks (7.16)

Your wheelchair is equipped with two sets of wheel locks, that operate from the left or the right side. They are applied directly against the tires. To engage, press the wheel-lock lever forward against the stop (fig. 78.2). To release the wheel, pull the lever back to the original position.

Braking power will decrease with:

- Worn tire tread
- Tire pressure that is too low
- Wet tires
- Improperly adjusted wheel locks.

The wheel locks have not been designed to be used as brakes for a moving wheelchair. The wheel locks should therefore never be used to brake a moving wheelchair. Always use the hand-rims for braking. Make sure that the interval between the tires and wheel locks complies with given specifications. To readjust, loosen screw and set the appropriate interval. Tighten screw (see the page on torque).

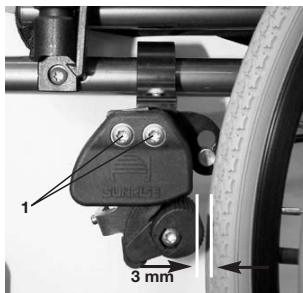78.1

78.2

Mounting the wheel lock too close toward the wheel, will result in a higher effort to operate. This might cause the brake lever extension to break! Leaning on to the brake lever extension while transferring will cause the lever to break! Splashing water from tires might cause the wheel lock to malfunction.

Options – Footplates

Footplates (7.21)

There are two types of footplates: platform and divided. They can be flipped up to facilitate getting in and getting out of your wheelchair.

Lower Leg Length

By removing the setscrews (1), the footrest can be adjusted to fit any lower leg length. Remove setscrews, set tubing with footrest in the desired position, and reinsert and tighten the setscrews (see the page on torque). Make sure that the composite plastic pad is in the proper position under the setscrew. A minimum interval of 2.5 centimeters from the ground should always be maintained.

13.1

Divided and Platform Footplates (7.22)

Footplates can be flipped up to facilitate getting in and getting out of your wheelchair.

They can also be tilted to six different angles relative to a level surface. Tighten screw (1) firmly on the outer side. By removing the clips (2), the footplate can be adjusted to three different positions toward both the front and rear. Loosen the adjustment screw (3) to change the horizontal position of the footplate. For this purpose, the footplate must be flipped up. When finished, make sure that all screws have been properly tightened (see the page on torque). A minimum interval of 2.5 centimeters from the ground should always be maintained.

15.1

15.2


15.3

Options – Footrests

Footrests and Latch (7.23)

Install your footrest by inserting it with the footplate(s) pointing toward the outside. Then rotate it toward the inside until it locks into place (1). To remove the footrest, depress the latch (2), swing the footrest toward the outside and remove it. Make sure that the footrest is properly locked into place.

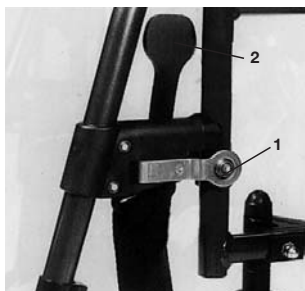16

CAUTION:

The footrests should not be used to lift or to carry the wheelchair with an occupant.

Little effort is needed to adjust the footrest by means of an eccentric hex nut. Adjustments can be changed as follows: Remove screw. Adjust the eccentric hex nut. Check if the locking mechanism engages and disengages properly. To make fast again, tighten the screw firmly while using a wrench to hold the eccentric hex nut securely (see the page on torque).

Flip-up footrest (standard or with length adjustment): (7.24)

To elevate:

Pull the footrest upwards and adjust it to the desired height. You will hear the footrest click into place.

To lower:

Lift your lower leg and press the release lever to unlock the footrest. You may now lower the footrest. As soon as you let go of the lever, the footrest will click back into place.

17

CAUTION:

Keep hands clear of the adjustment mechanism between the frame and the movable parts of the footrest while elevating or lowering the footrest.

CAUTION:

Footrests are not to be used for lifting or carrying the wheelchair with an occupant.

Options – Headrest

Headrest (7.35)

The headrest can be raised and moved both forwards and backwards. To do this, simply loosen the screw (1 or 2), adjust to the desired position, and tighten the screw (see the page on torque).

21

Options – Casters

Casters, Caster Plates, Forks (7.40)

The wheelchair may veer slightly to the right or left, or the casters may wobble. This may be caused by the following:

- Forward and/or reverse wheel motion has not been set properly.
- The camber has not been adjusted properly.
- Caster and/or rear wheel air pressure is incorrect; wheels do not turn smoothly.

22.1

The wheelchair will not move in a straight line if the casters have not been properly adjusted. Casters should always be adjusted by an authorized dealer. The caster plates must be readjusted, and the wheel locks must be checked any time the rear wheel position has been altered.

Easy^{max}

To ensure the forks are both set at the same angle, simply count the exposed splines on each side.

When the correct caster angle is found, the splines form a tight fit and provide 16° of adjustment in 2° increments.

Use this flat edge to square the fork.

The patented design allows the caster fork to rotate so it can remain square to the ground when the frame angle is changed.

191

192

193

194


195

Options – Axle Plates

Axle Plates (7.50)

The position of the body's center of gravity with respect to the rear wheels is critical for ensuring easy and pleasant use of the wheelchair. Several positions are possible: moving the axle plate (1) backwards along the holes drilled in the frame (2) will make the chair more stable, but it will not be as maneuverable as it would be if the axle plate were positioned farther to the front. You may adjust the seat height by moving the axle sleeve (3) in the axle plate. When finished, make sure that all screws have been properly tightened (see the page on torque).

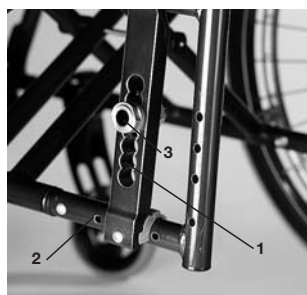27

CAUTION:

It may be necessary to readjust both the wheel locks and casters.


Options – Cambers

Cambers (7.51)

The camber may be set between 1° and 4°. This makes it possible to increase the chair's lateral stability as well as its maneuverability. Remove the axle plate screws (1). Adjust the adapter (2) to the desired angle. Tighten the screws. (see the page on torque). Raising the camber also raises the overall width of the wheelchair (by about 1 cm per degree).

28.1

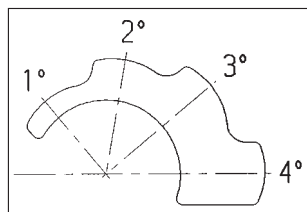28.2

CAUTION:

It may be necessary to readjust both the wheel locks and casters.

Options – Backrests

Angle-Adjustable Backrests (7.70)

The backrest can be adjusted to 5 different positions (3° forwards, 0°, 5°, 10°, and 15° backwards). A plate (2) is located on the sideguard, allowing you to set the angle by removing the screws (1) (see the page on torque).

35

CAUTION: The screws (1)

connect the backrest directly to the frame; screwing them in too tightly may cause them to break. Conversely, screws that are too loose can be easily lost. In either case, injuries could result when using the push handles to ascend stairs.

For your own safety, we recommend that adjustments only be made by authorized dealers.

Variable Angle-Adjustable Backrests (7.71)

You can adjust the angle of the backrest to meet your individual needs by pulling on the cord (1). The mechanism automatically snaps into place when you release the cord.

36

Adjustable Backrest Upholstery (7.72)

The sling can be adapted to any degree of tension by adjusting the Velcro fasteners. The padding in the backrest can be removed through an interior opening, allowing you to cushion the backrest according to your individual needs.

37


Options – Amputee Axle Plate

Amputee Axle Plate (7.55)

The larger the wheel base is, the greater is the stability of the wheelchair. The wheel base can be extended by screwing the axle sleeve (1) onto the frame at the desired height. Make sure that the steel pin (2) is mounted correctly in the hole drilled in the frame for this purpose. When finished, make sure that all screws have been properly tightened (see the page on torque).

32.1

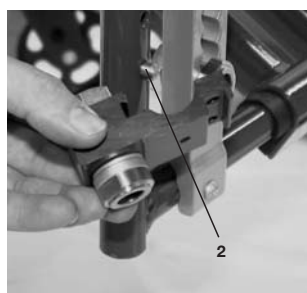32.2


CAUTION:

Wheel locks must be adjusted for the new position.

Options – Sideguards

Desk Sideguard, Flip-Up, Detachable with Short or Long Armrests, Height-Adjustable (7.81)

Adjust the height of the armrest as follows: Slide the knob (1) forward until it stops. Adjust the armrest (2) to the desired height. Slide the knob toward the back and push the armrest (2) down until you can hear it click into place. Always make sure that the sideguards have locked into place properly. To flip the armrest up, lift the handle (3) to unlock the sideguard.

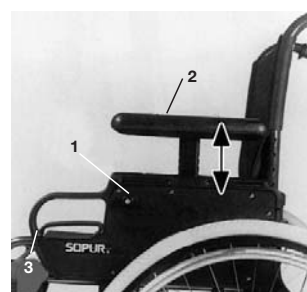45.1

CAUTION:

Neither the sideguards nor the armrests are to be used for lifting or carrying the wheelchair.

45.2

Options – Sideguards

Sideguards - Clothing Guards (7.82)

Clothing guards prevent water from splashing up and soiling clothing. Their position relative to the rear wheels can be set by adjusting the sideguards. To do this, remove the screws (1 and 2). After adjusting the sideguards to the desired position, replace and retighten the screws (see the page on torque).

46

CAUTION:

Neither the sideguards nor the armrests are to be used for lifting or carrying the wheelchair.

Sideguard, Aluminum, Bolt-On with Cold-Weather Protection (7.83)

Clothing guards (optional) prevent water from splashing up and soiling clothing.

Their position relative to the rear wheels can be set by adjusting the sideguards.

To do this, remove the screws (1 and 2). The armrests (optional) can be adjusted to the desired height after the screws (3) have been removed. When finished, make sure that all screws have been properly tightened (see the page on torque).

47


CAUTION:

Neither the sideguards nor the armrests are to be used for lifting or carrying the wheelchair.

Options – Push Handles

Height-Adjustable Push Handles (7.90)


These handles are secured with pins to prevent them from sliding out unintentionally. Opening the quick-release lever (1) makes it possible to adjust the push handles to meet your individual needs. As you move the lever, you will hear a locking mechanism; you may now easily position the push handle as desired. The nut (2) on the tension lever determines how tightly the push handles are clamped into place. If the nut is loose after adjusting the tension lever, the push handle will also be too loose. Turn the push handle from side to side before use to make sure that it is clamped securely enough into place. After adjusting handle height, always clamp the tension lever (1) securely into place. If the lever is not secure, injuries could result when lifting over obstacles.

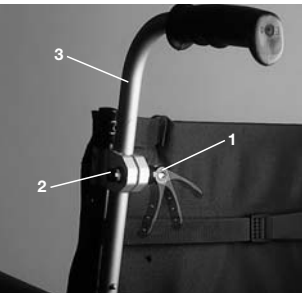49

Fold-Down Push Handles (7.91)

If the push handles are not in use, they can be folded down by depressing the button (2). When they are needed again, simply flip them back up until they click into place.

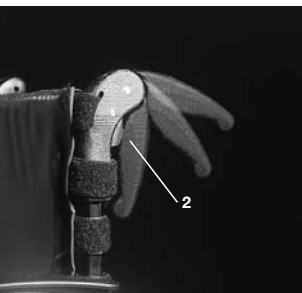50

Options – Lap belt instructions

Before using your wheelchair ensure the lap belt is worn.

The lap belt must be checked on a daily basis to ensure it is free from any obstruction or adverse wear.

Always make sure that the lap belt is correctly secured and adjusted prior to use. Too loose a strap could cause the user to slip down and risk suffocation or cause serious injury.

The lap belt is fitted to the wheelchair as shown in the illustrations. The seat belt comprises 2 halves. They are fitted using the existing seat stay retaining bolt fitted through the eyelet on the belt. The belt is routed under the rear of the side panel. (Fig. 24)

Adjust the belt position so buckles are in the centre of the seat. (Fig. 25)

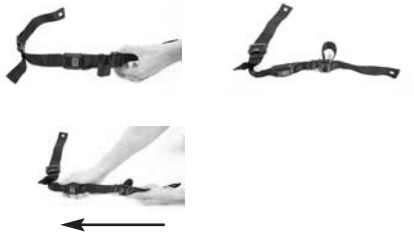
Fig. 24

Fig. 25

Adjust lap belt to suit the user's needs as follows:

To reduce the belt length	To increase the belt length
	
Feed free belt back through male buckle and slide adjusters. Ensure belt is not looped at male buckle.	Feed free belt through slide adjusters and male buckle to provide more belt length.

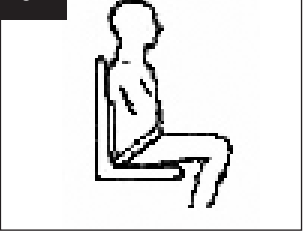
When fastened check space between the lap belt and user. When correctly adjusted it should not be possible to insert more than the fat of the hand between the lap belt and the user. (Fig. 26)

Fig. 26

The lap belt should be fixed so that the belt sits at an angle of 45 degrees across the users pelvis. The user should be upright and be as far back as possible in the seat when correctly adjusted. The lap belt should not allow the user to slip down in the seat. (Fig. 27)

Fig. 27

	
To fasten buckle: Firmly push male buckle into female buckle.	To release belt: Press exposed sides of male buckle and push towards centre whilst gently pulling apart.

 If in doubt about the use and operation of the lap belt then ask your healthcare professional, wheelchair dealer, carer or attendant for assistance.

Advice to client

The lap belt must only be fitted by an approved Sunrise Medical dealer / agent. The lap belt should only be adjusted by a professional, or a Sunrise Medical approved dealer / agent. The lap belt must be checked on a daily basis to ensure it is adjusted correctly and free from any obstruction or adverse wear. Sunrise Medical does not encourage the transportation of any person in a vehicle using this lap belt as a method of restraint. Please see Sunrise Medical transit booklet for further advice on transportation.

Maintenance:

Check lap belt, and securing components, at regular intervals for any sign of frays or damage. Replace if necessary.

NOTE:

The lap belt should be adjusted to suit the end user as detailed above. Sunrise Medical recommend that the length and fit of the belt be checked on a regular basis to reduce the risk of the end user inadvertently re-adjusting the belt to an excessive length.

Options – Anti-Tip Tubes

Anti-Tip Tubes (7.93)

Anti-tip tubes (1) provide additional safety for inexperienced users when they are still learning how to operate their wheelchairs. Anti-tip tubes (1) prevent a wheelchair from tipping over backwards. Pushing the anti-tip tubes (1) will move them downwards; they can also be swung forward. An interval of 3 cm to 5 cm should be maintained between them and the ground. You must swing the anti-tip tubes forward when going up and down large obstacles (such as a curb) to prevent them from touching the ground.


53

Options – Standard Sling

Standard Sling (7.101)

The sling is equipped with a Velcro fastener on one side; this permits stepless adjustment of the sling. To readjust the sling, first fold up the wheelchair slightly. Then remove the screws (1) and pull the front plugs (2) forward out of the frame. The sling (3) can then be slid off of the frame. By unfastening the Velcro, the sling can be readjusted. Reverse this procedure to reinstall the sling. When finished, make sure that all screws have been properly tightened (see the page on torque).

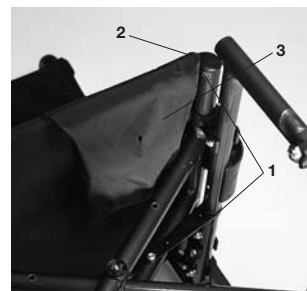54

CAUTION:

For proper safety, at least 50 % of the Velcro surfaces must make contact with each other at all times.

Options – Crutch Holder

Crutch Holder (7.110)

This device permits crutches to be transported directly on a wheelchair. It has a Velcro loop (1) to fasten crutches or other aids.

CAUTION:

Never try to use or even remove the crutches or other aids while moving.

58

Options – Tray

Tray (7.111)

The tray provides a flat surface for most activities. Before using a tray, it first has to be adjusted once to the width of the seat by an authorized dealer. The user must be sitting in the wheelchair during the trial run.

60

Options – Stabilizing Bar

Stabilizing Bar (7.112)

This bar is used to stabilize the backrest. To be able to fold up your wheelchair, unscrew the knobs (1) from the push handles and remove the stabilizing bar; now you can easily fold up the wheelchair. To reinstall the stabilizing bar, the knobs must be reinserted and tightly screwed so that the bar sits firmly in position.

61

Options – Transit Wheels

Transit Wheels (7.113)

Transit wheels should be used wherever your wheelchair would be too wide if the rear wheels were used (e.g., in airplanes, buses, etc.). After the rear wheels have been removed with the aid of the quick-release axles, the transit wheels can immediately be used to continue riding. The transit wheels are mounted so that they are approx. 3 centimeters above the ground when not in use. They are thus out of the way when riding, transporting, or when tipping to pass over obstacles (e.g., curbs, steps, etc.).

CAUTION:

Your wheelchair does not have any wheel locks when the transit wheels are being used.

49.1

49.2

Tires and Mounting

Tires and Mounting (8.1)

Always make sure you that you maintain the correct tire pressure, as this can have an effect on wheelchair performance. If the tire pressure is too low, rolling resistance will increase, requiring more effort to move the chair forward; low tire pressure also has a negative impact on maneuverability. If the tire pressure is too high, the tire could burst. The correct pressure for a given tire is printed on the surface of the tire itself. Tires can be mounted the same way as an ordinary bicycle tire. Before installing a new inner tube, you should always make sure that the base of the rim and the interior of the tire are free of foreign objects. Check the pressure after mounting or repairing a tire. It is critical to your safety and to the wheelchair's performance that regulation air pressure be maintained and that tires be in good condition.

Possible Problems

Wheelchair pulls to one side (9.1)

- Check tire pressure
- Check to make sure wheel turns easily (bearings, axle)
- Check angles of casters
- Check to make sure both casters are making proper contact with the ground

Casters begin to wobble (9.2)

- Check angle of casters
- Check to make sure all bolts are secure; tighten if necessary (See the page on torque)
- Check to make sure both casters are making proper contact with the ground

Wheelchair / Cross-tube assembly does not snap into position in the seat saddle (9.3)

- Chair is still new, i.e., the seat or backrest upholstery is still very stiff. This will improve with continued use.

Wheelchair is difficult to fold up (9.5)

- Adjustable backrest upholstery is too stiff; loosen it accordingly.

Wheelchair squeaks and rattles (9.6)

- Check to make sure bolts are secure; tighten if necessary (see the page on torque).
- Apply small amount of lubrication to spots where movable parts come in contact with one another

Wheelchair begins to wobble (9.7)

- Check angle at which casters are set
- Check tire pressure
- Check to see if rear wheels are adjusted differently

Maintenance

Maintenance (9.9)

- Check the tire pressure every 4 weeks. Check all of the tires for wear or damage.
- Check the brakes approximately every 4 weeks to make sure that they are working properly and easy to use.
- Change tires as you would an ordinary bicycle tire.
- All of the joints that are critical to safe using your wheelchair safely are self-locking nuts. Please check every three months to make sure that all bolts are secure (See the page on torque). Safety nuts should only be used once and should be replaced after multiple uses.
- Please use only mild household cleansers when your wheelchair is dirty. Use only soap and water when cleaning the seat upholstery.
- If your wheelchair should ever get wet, please dry it after use.
- A small amount of sewing-machine oil should be applied to quick-release axles approximately every 8 weeks. Depending on the frequency and type of use, we recommend taking your wheelchair to your authorized dealer every 6 months to have it inspected by trained personnel.

CAUTION:

Sand and sea water (or salt in the winter) can damage the bearings of the front and rear wheels. Clean the wheelchair thoroughly after exposure.

Nameplates / Guarantee

Nameplates (11.1)

The nameplate is located on either the cross-tube assembly or the transverse frame tube, as well as on the back page of the owner's manual. The nameplate indicates the exact model designation and other technical specifications. Please provide the following pieces of information whenever you have to order replacement parts or to file a claim:

- Serial number
- Order number
- Month / Year

Guarantee

You have purchased a high-quality Sopur product. As a sign of our gratitude, we are providing you with a 5-year guarantee on all frame components and for the cross-tube assembly. We are not responsible for any damage resulting from inappropriate or unprofessional installation and/or repairs, neglect, wear, or from changes in any wheelchair components caused either by the user or by third parties. In such cases, this guarantee shall be considered null and void.

Custom wheelchairs cannot be exchanged.

Technical Data – Possible seat-height adjustments

Possible seat-height adjustments

WARNING: Measured without seat cushions

Type of caster	Type of fork	Front seatheight in cm	Rear seat-height in cm				
			22" rear wheel only		22" and 24" rear wheel		24" rear wheel only
4" solid rubber	Short (4-hole)	43	40	41	43		
	Short (4-hole)	45	40	41	43	44	
	Short (4-hole)	46	40	41	43	44	46
	Long (5-hole)	47	40	41	43	44	46 47
	Long (5-hole)	48	40	41	43	44	46 47 48
5" solid rubber	Short (4-hole)	44	40	41	43	44	
	Short (4-hole)	45	40	41	43	44	
	Short (4-hole)	46	40	41	43	44	46
	Long (5-hole)	47	40	41	43	44	46 47
	Long (5-hole)	48	40	41	43	44	46 47 48
	Long (5-hole)	49	40	41	43	44	46 47 48 49
6" solid rubber	Long (5-hole)	47	40	41	43	44	46 47
	Long (5-hole)	48	40	41	43	44	46 47 48
	Long (5-hole)	49	40	41	43	44	46 47 48 49
	Long (5-hole)	50	–	41	43	44	46 47 48 49
6" pneumatic	Long (5-hole)	47	40	41	43	44	46 47
	Long (5-hole)	48	40	41	43	44	46 47 48
	Long (5-hole)	49	40	41	43	44	46 47 48 49
	Long (5-hole)	51	–	–	43	44	46 47 48 49 51
6" soft wheel	Long (5-hole)	46	40	41	43	44	46
	Long (5-hole)	48	40	41	43	44	46 47 48
	Long (5-hole)	49	40	41	43	44	46 47 48 49
	Long (5-hole)	50	–	41	43	44	46 47 48 49
7" pneumatic/solid rubber	Long (5-hole)	48	40	41	43	44	46 47 48
	Long (5-hole)	49	40	41	43	44	46 47 48 49
	Long (5-hole)	50	–	41	43	44	46 47 48 49
	Long (5-hole)	52	–	–	43	44	46 47 48 49 51

Possible seat-height adjustments Easy^{max}

WARNING: Measured without seat cushions

Type of caster	Seat height front in cm	Seat height rear in cm
6"	50	44, 45, 46, 47
	51	44, 45, 47, 48, 50
	52	45, 47, 48, 50
	53	47, 48, 50
5"	48	43, 45, 46, 47
	49	43, 44, 46, 47
	50	43, 44, 46, 47, 48, 49
	51	44, 45, 46, 47, 48, 49, 50
	52	47, 48, 50
4"	46	43, 45
	47	43, 45, 46
	48	43, 44, 45, 46, 48
	49	43, 44, 45, 46, 47
	50	43, 44, 46, 48, 49
	51	47, 48, 50

Technical Data

Technical Data (Easy 200-300)

Total width:

- with standard wheels, incl. Hand rims, close mounted
- in combination with desk sideguards: SW + 20 cm
- in combination with aluminium sideguards: SW + 18 cm
- in combination with clothing protection, composite: SB + 19 cm
- A super-narrow handrim attachment reduces the overall width by 2 cm
- Wheels with drum brakes widen the chair by 2 cm

Dimensions, folded:



- without rear wheels: approx. 30 cm
- with rear wheels: approx. 25 cm

Weight in kilograms:

- Footrests Easy 300: 1,35 kg
- Footrests Easy 200: 0,93 kg
- Chair with no sideguards, footrests or rear wheels: 7,60 kg
- Desk-Sideguards: 3,10 kg


Maximum weight limit:

- Approved for user weights up to 120 kg/265 lbs.

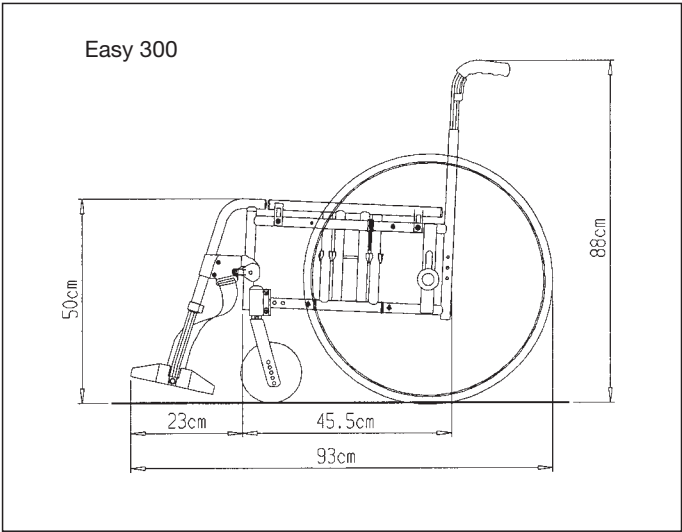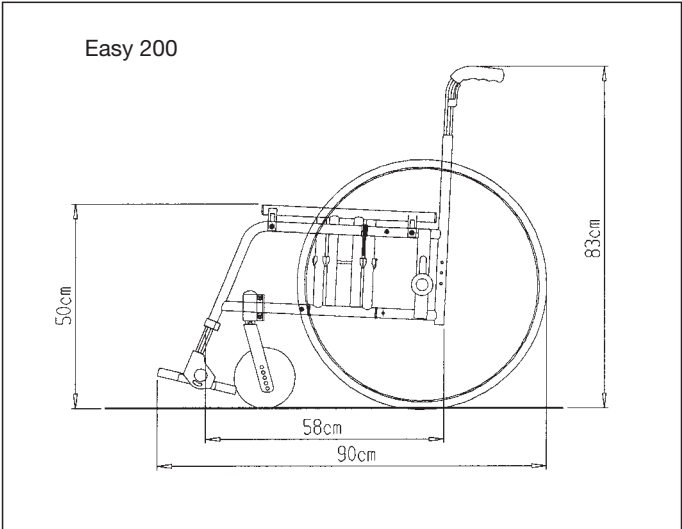
mount position axle plate	mount position brake bracket	Brake bracket type
short frame right		
1	1	left
2	2	left
3	2	left
4	1	left
5	1	left
long frame right		
1	2	right
2	2	right
3	1	right
4	1	right
5	2	left
short frame left		
1	1	left
2	2	right
3	2	right
4	1	right
5	1	right
long frame left		
1	2	left
2	2	left
3	1	left
4	1	left
5	2	right


 = presetted

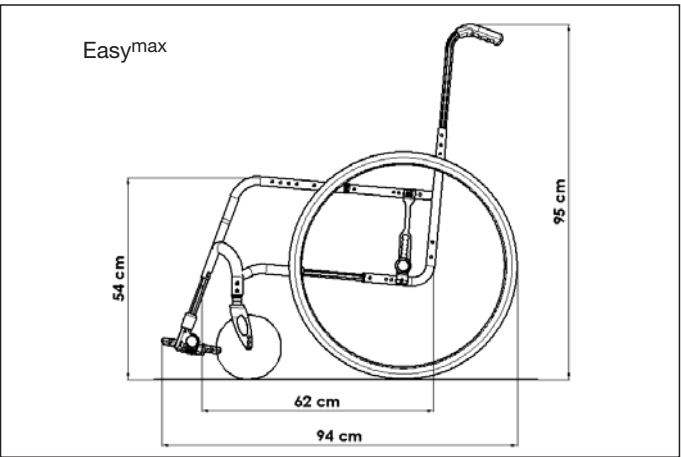
199

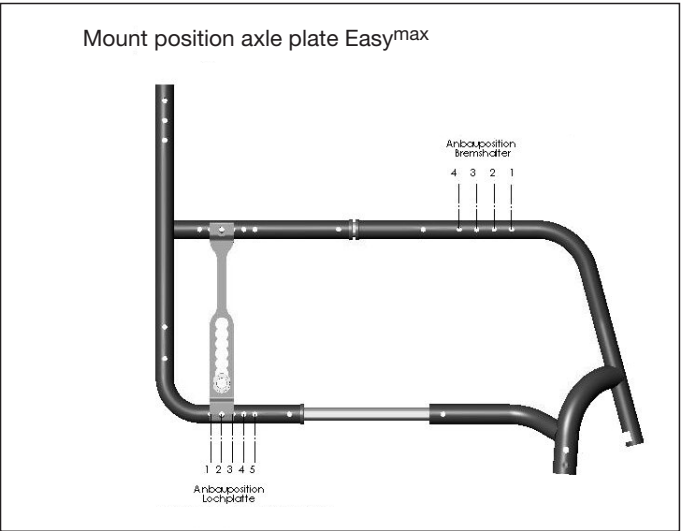196

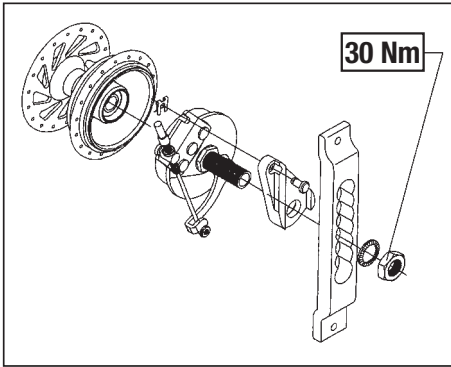
64

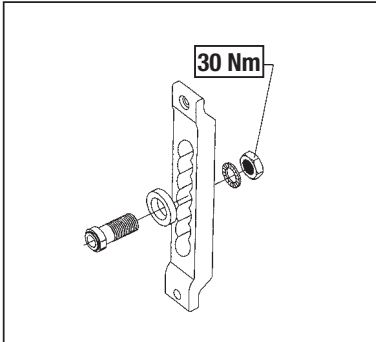
69

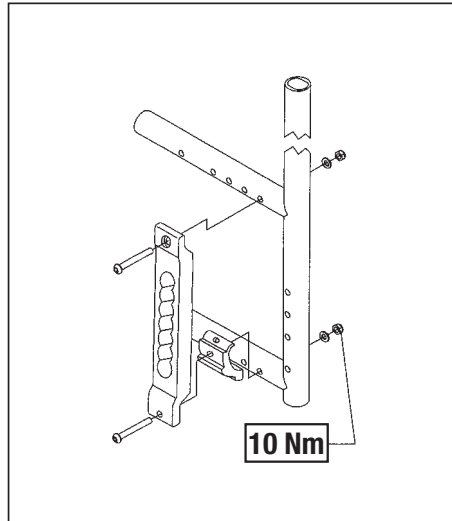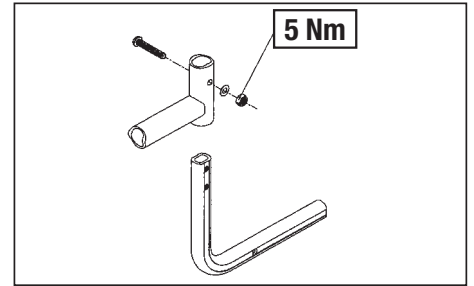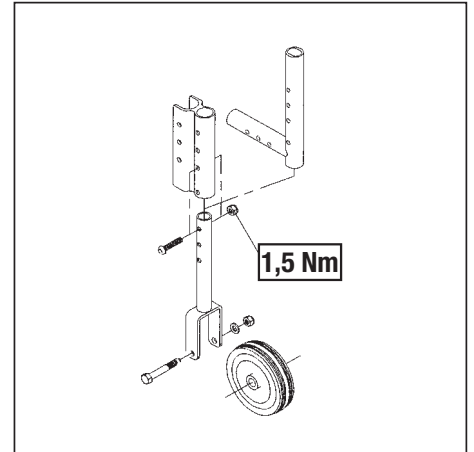
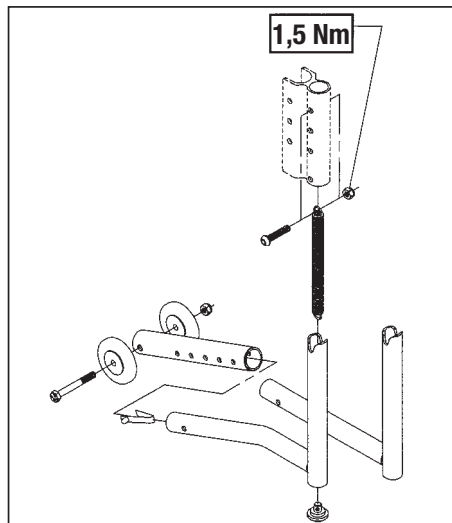
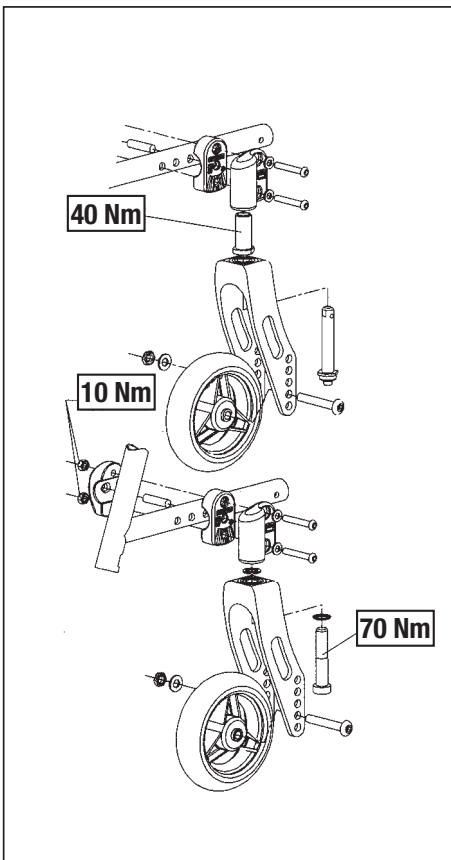
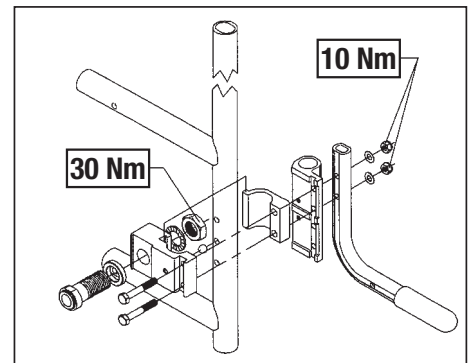
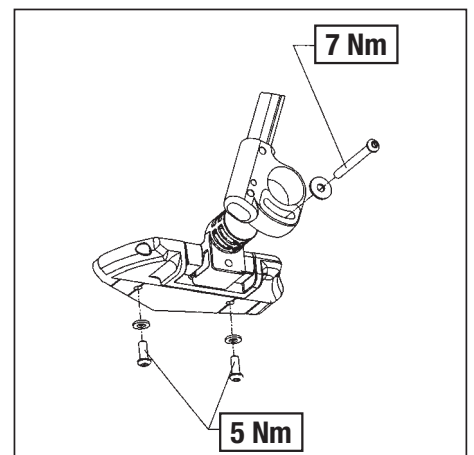
The torque for the M6 screw is 7 Nm, unless otherwise specified.

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Sunrise Medical AB
Sunrise Medical B.V.
G2 A/S